

Transform  
your body  
with the power  
of bioquantic  
synergies.



mya 6.0  
Bioquantic synergies

ADVANCED TECHNOLOGY,  
UNLIMITED POSSIBILITIES.



TECHNOLOGY  
OF EXCELLENCE  
MADE IN ITALY



 **mya 6.0**  
CESARE  
QUARANTA Bioquantic synergies

## 6 technologies for results that are 6 times more effective

**MYA 6.0** is the latest **hand-free total body** device from Cesare Quaranta, the only technology on the market offering **truly outstanding results**. By harnessing the **synergy between six technologies**, customers' beauty and wellbeing goals are naturally within reach thanks to **rapid and effective bodily and facial treatments**.

The cellular and tissue bio-vitalisation and bio-regeneration processes are achieved thanks to the LLLT laser which enhances the effects of the other technologies. MYA 6.0 brings together the 6 best beauty care technologies in a single device

- LLLT LASER
- ELECTROSTIMULATION
- ULTRASOUND
- ELECTROPORATION
- PHOTOBIMODULATION
- BIA TEST

# What makes our products different

## EXPERIENCE THE NEW MYA 6.0 BIOQUANTIC SYNERGY METHOD

Building on our experience in the development and production of our patented ultrasound and electrostimulation pads, we have incorporated LLLT laser technology into the pads to enhance and maximise benefits and results. LLLT LASER - LOW LEVEL LASER THERAPY technology penetrates deeply, directly to the hypodermis and adipose tissue, at a wavelength of 650nm or 780nm.

The bioquantum energy sets off a domino effect of reactions which intensify and boost the effects of the energy applied to the tissues, effectively countering imperfections for results better than ever seen before.

## ULTRASOUND AND LLLT LASER THE PERFECT SYNERGY FOR 6 TIMES MORE EFFECTIVE RESULTS

Ultrasound and LLLT laser work in synergy to enhance the effects of electrostimulation on the target tissue, depending on the chosen program. Results are enhanced by up to 6 times in a single treatment. The powerful action of LLLT laser technology on adipose tissue accelerates weight loss by ensuring the enhancement of tissue tone through the stimulation of collagen.

## OUR FULLY REDESIGNED PADS

**OVAL-SHAPED AND ERGONOMIC:** the new ultra-ergonomic shape is perfect for shock absorption.

**LARGER:** to maximise the contact and treatment area.

**NEW COLOURS:** white and blue for an attractive, innovative appeal.

**OPTIMISED TO REFLECT THE LLLT LASER BEAM:** ensuring the complete safety of the operator and customer.

+ LLLT LASER

+ ELECTROSTIMULATION

+ HIGH FREQUENCY ULTRASOUND





miva 60



4

4

3

3

2

2

1


1

## 6.0 BENEFITS FOR BODY SCULPTING

1. SCULPTING
2. SLIMMING
3. LYMPHATIC DRAINAGE
4. MUSCLE TONING
5. TISSUE TONING
6. CELLULITE TREATMENT

While the same pad also treats:

- LOCALISED ADIPOSE DEPOSITS;
- CELLULITE;
- POOR MUSCLE TONE;
- LYMPH STASIS;
- OVERALL WELLBEING

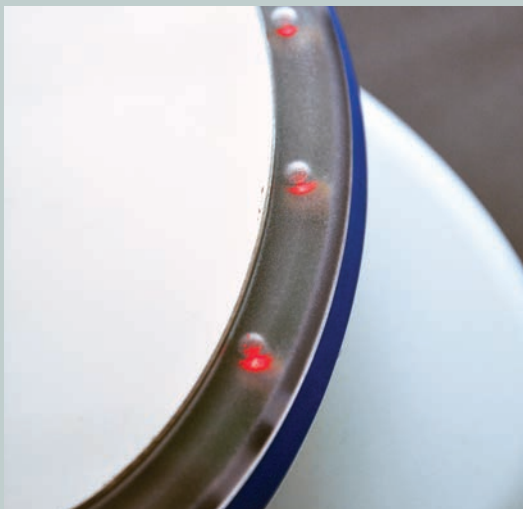

### CUSTOMISED 15-TREATMENT REMODELING PROGRAM.

The customer **lost 6 cm off the abdominal area**, achieving a slimmer and more toned shape.



### CUSTOMISED 12-TREATMENT REMODELING PROGRAM.

Visible reduction in oedematous and fibrous thigh cellulite.



# A one-of-a-kind new method for unparalleled results

## CASE STUDY 1

### BEFORE TREATMENT:

- + fibrous cellulite
- + localised adipose deposits on the hips, buttocks, the backs of the thighs and the upper outer thighs
- + small amount of soft cellulite on the inner thigh

### ANTHROPOMETRIC MEASUREMENTS

12 bi-weekly treatments.

Total treatment time: approximately 2 months.

|       | Pre-treatment | After 6 sessions | After 12 sessions |
|-------|---------------|------------------|-------------------|
| HIPS  | 103 cm        | 100 cm           | 96 cm             |
| THIGH | 65 cm         | 63 cm            | 60 cm             |

**WEIGHT: 58 Kg (UNCHANGED)**

Adipose tissue replaced with lean tissue.

### AFTER TREATMENT:

- + the orange peel effect due to fibrous cellulite is almost entirely gone
- + visible loss of centimetres off the hips and upper outer thighs
- + significant improvement in muscle tone
- + the soft cellulite on the inner thigh is gone



## CASE STUDY 2

### BEFORE TREATMENT

- + significant abdominal adipose deposits
- + oedematous and fibrous cellulite

### ANTHROPOMETRIC MEASUREMENTS

15 bi-weekly treatments

WEIGHT DECREASED FROM 63 Kg TO 60 Kg

|         | Pre-treatment | After 15 sessions |
|---------|---------------|-------------------|
| WAIST   | 97 cm         | 88 cm             |
| ABDOMEN | 112 cm        | 103 cm            |
| THIGH   | 60 cm         | 57 cm             |

### AFTER TREATMENT

- + significant decrease in the abdominal adipose deposits
- + visible reduction in oedematous and fibrous cellulite



# Bioquantic synergies 6.0 technologies

6 TECHNOLOGIES WORKING IN SYNERGY FOR BIOQUANTIC BENEFITS OFFERING RESULTS UNLIKE ANYTHING SEEN BEFORE

## 01. LLLT LASER

The source of the MYA 6.0 bioquantic synergies: each pad has 6 LLLT lasers with wavelength of 635nm of 672mW operating on an intermittent duty cycle and the activation and relaxation phases are alternated to maximise the effects. LLLT laser technology penetrates the skin barrier and triggers photochemical responses in the treated tissues, inducing metabolic reactions in the adipocytes and triggering a cascade of reactions.

## 02. ELECTROSTIMULATION

Electrostimulation uses deep impulses on the face and body which lead to muscle and skin lifting, improving oxygenation and tissue revitalisation and facilitating metabolic exchange and the elimination of toxins. The controlled emission of electrical currents with different waveforms at frequencies ranging from 33 Hz to 5,000 Hz enables effective beauty treatments which can be tailored to the customer's individual needs.

## 03. ULTRASOUND

Ultrasound acts on cellulite by combining thermal and localised micromechanical effects. In addition to massage and localised heat production, it offers improvements in vascularisation and molecular dissociation.

## 04. BIA TEST

Bioelectrical impedance analysis assesses a customer's body composition. Its sinusoidal current with three-phase frequency enables accurate, reliable measurements of the percentage of lean tissue, body fat and overall intracellular and extracellular liquids. This technology analyses a customer's condition before treatment and their improvement over the course of the treatment.

## 05. ELECTROPORATION



Electroporation works by evenly distributing the cosmetic product over the skin and facilitating the complete absorption of its active ingredients.

The electroporation handpiece improves tissue permeability by making small, reversible structural modifications to the stratum corneum

## 06. PHOTOBIO-MODULATION

Photobiomodulation can be seamlessly incorporated into wellbeing treatments, complementing traditional beauty treatments for an enhanced overall result. Treatment involves the use of different coloured lights to improve specific areas of the body: imperfections are considered inharmonious vibrations of the cells which can be rebalanced by channelling the energy - and thus the vibration - of one or more colours.



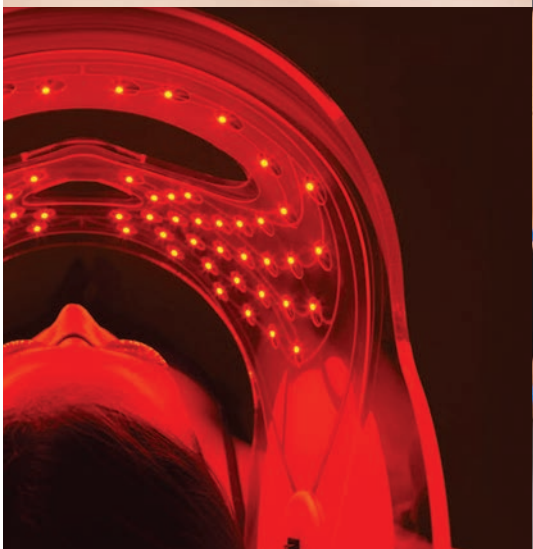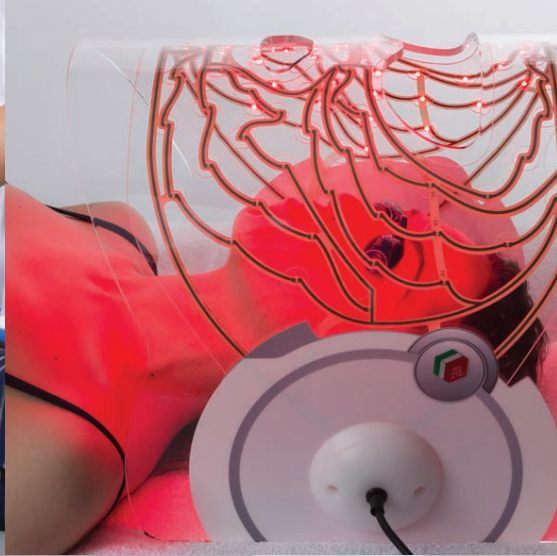


# 31 programs, endless possibilities.

**Research & Development is a foundational value of our company.**

Our team of experts carefully designs and develops programs and treatments for our devices with the needs of beauty care professionals and end customers in mind. The MYA 6.0 programs have been developed following an in-depth analysis of bioquantic synergies, and the LLLT laser - the source of these synergies - **has been incorporated into all programs to boost the effectiveness of the pads and to enhance results.** The range of programs available has also been expanded and diversified over the years, offering a comprehensive range of treatments:

- 2 programs for electrostimulation with electrodes for the facial area
- 11 programs for electrostimulation with electrodes for the body
- 13 programs for ultrasound / LLLT laser + electrostimulation within the same body pad
- 1 program for the absorption of cosmetics and photobiomodulation
- 1 program for BIA test
- 1 chip card program
- 1 counter program
- 1 adjustments program

All programs can be tailored to the customer and the required treatment.

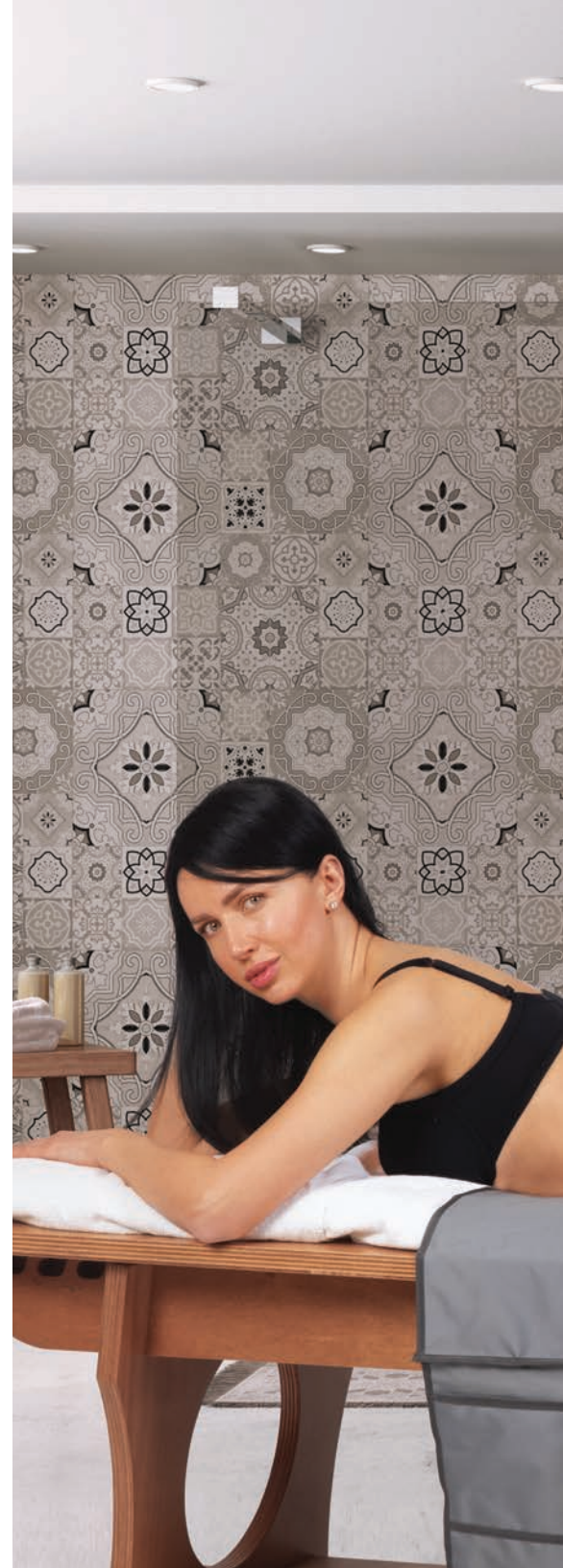



## TREAT YOUR FACE AND BODY AT THE SAME TIME WITH MYA 6.0

The device comes with an **ANTI-AGEING LED MASK**. The **photobiomodulation mask is equipped with 70 red LED lights** to counteract the signs of ageing and skin laxity. The stimulation is athermic, painless and can be used on all skin phototypes.

- Increased cell turnover at the epidermal level;
- Increased collagen production by fibroblasts;  
Improved skin texture;
- Decrease in fine lines;
- Skin tightening.



## A REVOLUTION IN THE FIELD OF PHOTOREJUVENATION: LED MASK + POLYNUCLEOTIDE SERUM

Facial skin laxity, dullness, lines and hyperpigmentation: treat face, neck and décolletage blemishes with a program designed to restore the skin's vitality and beauty. **The LED MASK combined with the effectiveness of a professional polynucleotide cosmetic mask is a revolutionary one-step photorejuvenation solution.** Ideal both for women and men, this combination boosts the cellular reaction and triggers intense tissue stimulation, offering exceptional results.



mya60



4

4

3

3

2

# Certified electro-medical quality



Cesare Quaranta is an Italian company that designs and produces electro-medical technologies.

The company has two ISO certifications:

## ISO 9001:2015:

Design, manufacture, service and placing on the market of electromedical devices for dermatology and aesthetic medicine. Design, manufacture, service and placing on the market of electronic equipment for the aesthetic sector. Marketing and sale of cosmetic products.

## ISO 13485:2016:

Design, manufacture service and placing on the market of electromedical devices for dermatology and aesthetic medicine.

All devices are produced in our laboratories in Chieri (Turin) using state-of-the-art tools and undergo rigorous testing.



## TECHNICAL DATA

### ELECTROSTIMULATION

|                       |  |
|-----------------------|--|
| Output stage          | square wave - Kotz currents                          |
| Square wave frequency | 1kHz modulated to 33/66/133/200 Hz                   |
| Kotz wave frequency   | 2500/4000/5000 Hz modulated to 4/6/8/12/50/75/100 Hz |
| Output current (max)  | 80 mA with load from 500 ohm                         |

### ULTRASOUND

|                      |                         |
|----------------------|-------------------------|
| Ultrasound frequency | 3 MHz                   |
| Densità di potenza   | < 2.8 W/cm <sup>2</sup> |

### LLLT LASER

|               |                 |
|---------------|-----------------|
| Wavelength    | 650 nm o 780 nm |
| Maximum power | 672 mW          |

### COSMETIC ABSORPTION AND PHOTOBIMODULATION

|                       |                                   |
|-----------------------|-----------------------------------|
| Output stage          | square wave                       |
| Square wave frequency | 1,6 kHz modulated to 300 mHz      |
| Output current (rms)  | 80 mA with load from 500 ohm      |
| Colours               | blue, white, yellow-orange, green |

### BIOIMPEDENTIOMETRY

|                           |  |
|---------------------------|--|
| Current for BIA           | 0,8 mA a 5 kHz, 50 kHz, 100 kHz  |
| Number of output channels | 8 right and 8 left channels for ultrasound, 4 right and 4 left channels for electrostimulation, individually adjustable for each output with analog and digital information on the display |
| Selection of treatments   | 13 preset programs for traditional electrostimulation<br>13 preset programs for Ultrasound, LLLT Laser and/or Electrostimulation.  |
| Total treatment time      | Preset or can be programmed from 5 to 60 minutes   |

|                   |                                  |
|-------------------|----------------------------------|
| Power supply      | 220/250 V. single-phase 50/60    |
| Absorbed power    | 190 VA MAX                       |
| Mains fuses       | 2 x 2 A type T 5 x 20            |
| Size              | 420 x 1330 x 470 mm              |
| Weight            | 22 kg                            |
| Temperature range | from 5° C to 40° C (when in use) |
| Product           | Made in Italy                    |

**NB:** The manufacturer reserves the right to make technical and design changes to improve the product without prior notice.

NB: The manufacturer reserves the right to make technical and design changes to improve the product without prior notice.

  
TECHNOLOGY  
OF EXCELLENCE  
MADE IN ITALY



C E S A R E  
Q U A R A N T A

**CESARE QUARANTA s.r.l - Unipersonale**  
Via Legnanino, 5 - 10023 Chieri (TO) Italy - +39 0119472132  
[www.cesarequaranta.it](http://www.cesarequaranta.it) - e-mail: [info@cesarequaranta.it](mailto:info@cesarequaranta.it)